

## Understanding changes in porosity as a function of overburden pressure in shales

Application Note 5



For the profitable development of shale reservoirs, it is critical to understand the natural fractures in the rock and how they may or may not interact with hydraulic fracturing. When performing core analysis, it is also important to assess whether the measured properties have been altered by fractures induced during coring.

A key indication of the quantity of the fractures can be obtained by measuring the amount of porosity contained in the fracture network. In shale formations there usually exists at least two distinct pore networks - one associated with the intergranular pore network and one associated with the pores within the organic rich kerogen. The fracture network, either natural or induced, could be thought of as a third pore network which at ambient pressure would be the largest in size. The  $T_2$  distribution is still able to distinguish the different pore networks because the organic pores are much smaller than their non-organic counterparts. In addition, the  $T_2$  distributions of the three pore networks should change differently as a function of confining pressure allowing further distinction of the three pore networks to be made.

In this application note, we describe techniques using NMR that can not only obtain the total porosity of shale samples but can also quantify the amount of porosity arising from the fractures. NMR measurements of the  $T_2$  relaxation time were performed at different confining pressures to quantify the porosity loss as confining stress increases.

### Method

The total porosity and  $T_2$  distribution as a function of confining pressure of four shale core samples was investigated using NMR. Each dry shale sample was saturated with brine via pressure saturation at 10,000 psi using a Phoenix Instruments pressure vessel. Following this saturation procedure, each shale sample was confined hydrostatically by fluorinert in an Oxford Instruments P5 overburden NMR probe. Once the pressure was stabilized, the cell was inserted into an Oxford Instruments **GeoSpec** 2/53 rock core analyzer. The porosity of each shale core sample was then tested as a function of confining pressure at 0, 1000, 2000, 3000, 4000 and 5000 psi. A  $T_2$  NMR acquisition scan of each sample, at each pressure, was used to measure the  $T_2$  pore distribution as a function of confining pressure. The  $T_2$  distribution of a fluid saturated rock is proportional to its pore size distribution. Fluid in smaller pores has a shorter  $T_2$  value while fluid in larger pores has a longer  $T_2$  value. The pore volume at each confining pressure was determined by summing the area under the  $T_2$  distribution curve. The NMR instrument had been calibrated with a sample of known volume so the NMR signal amplitude indicated the volume of fluid within the shale samples. Acquisition and analysis of the  $T_2$  data was achieved via Green Imaging Technologies' software.



*The Business of Science®*

# Understanding changes in porosity as a function of overburden pressure in shales



## Results

The pore-size distribution for each shale core sample tested at 5,000 psi of confining pressure can be seen in Figure 1. Shales 1 and 3 have a single sized pore network while shales 2 and 4 have three pore networks (labeled: micro, meso and macro) of varying size. For shale samples 2 and 4, the cumulative porosity of each network was estimated by summing the incremental porosity in certain  $T_2$  ranges corresponding to the different pore networks.

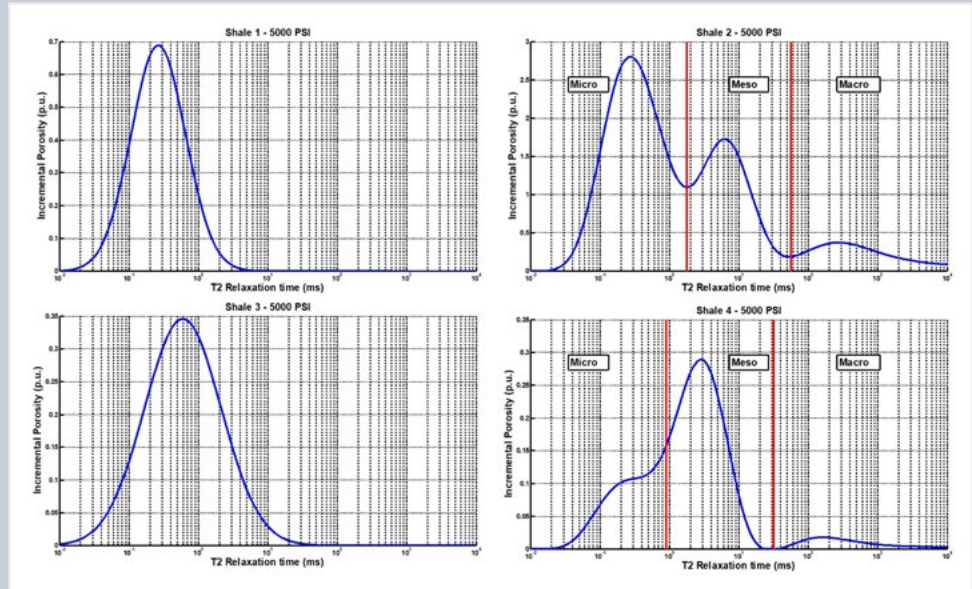


Figure 1: Pore distribution for shale samples tested at 5000 psi confining pressure

The left panel of Figure 2 shows the porosity of each pore network of shale 2 plotted as a function of confining pressure. This plot shows that the micro pore network makes up the majority of the total porosity of shale 2 (~2.5 p.u.) while the macro porosity (~0.5 p.u.) accounts for the least portion of the total porosity of shale 2. The center panel of Figure 2 shows the porosity normalized to 0 psi and plotted as a function of pressure.



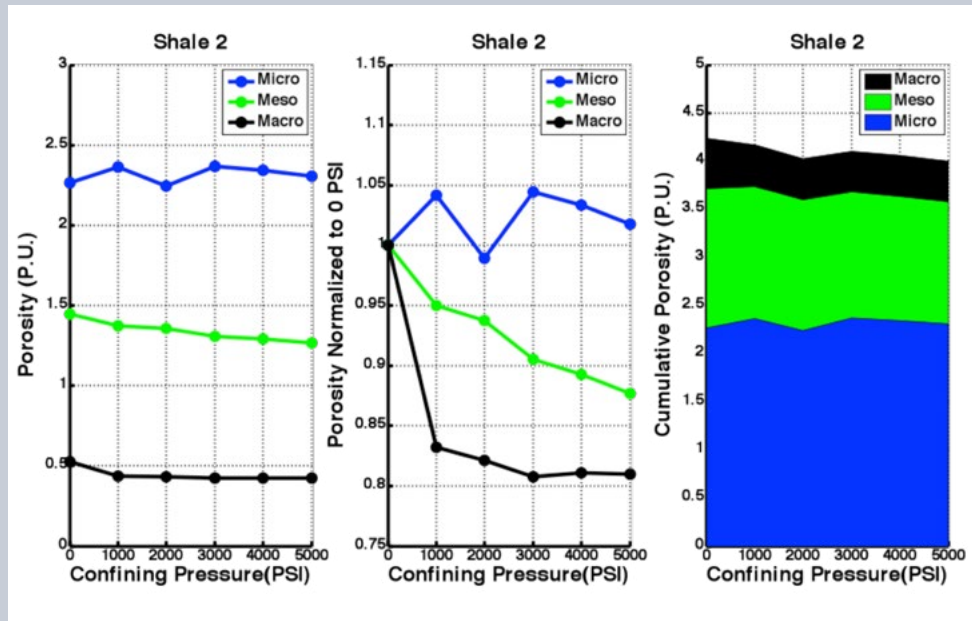



Figure 2: The porosity of the pore networks of shale 2 as a function of pressure

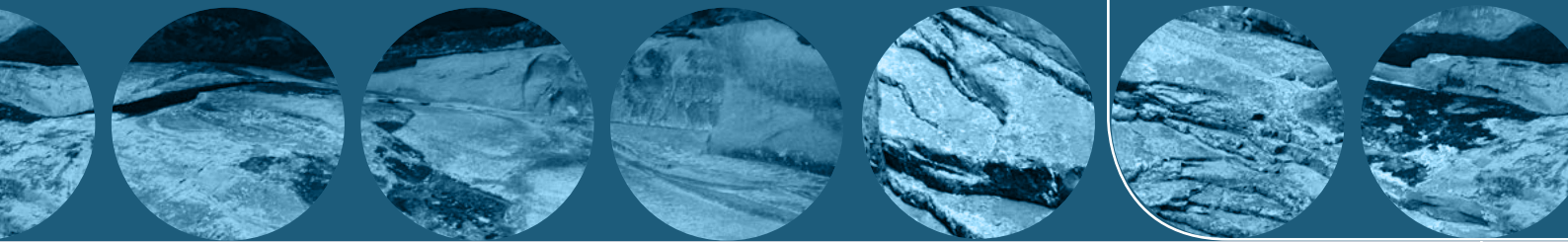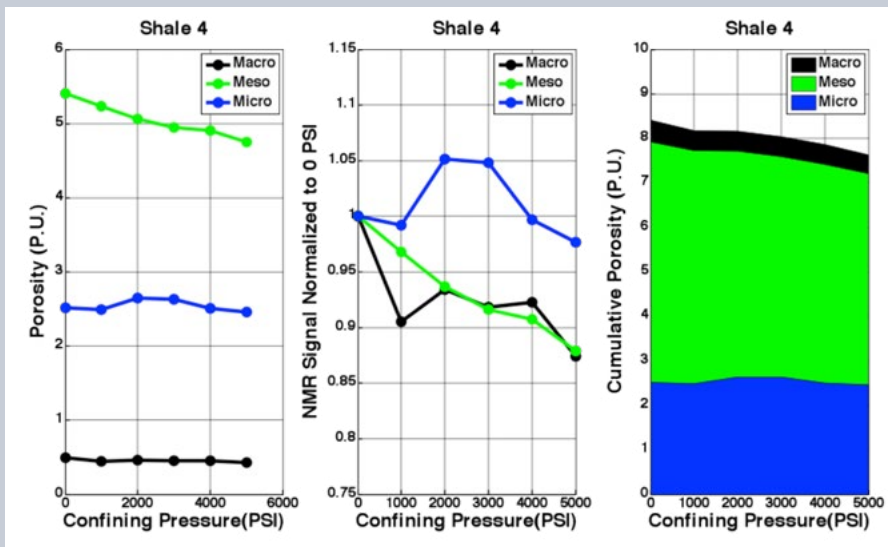

From this plot it is immediately obvious that each pore network behaves differently with confining pressure. The micro pore network (blue line) shows no compression with confining pressure. The meso pore network (green line) shows an approximately linear compression with confining pressure. The macro pore network shows a sharp decrease in porosity from 0 to 1,000 psi and then no further compression from 1,000 to 5,000 psi. This is consistent with the observations of Chhatre et al. [1]. In their study, there was a significant decrease in permeability with increasing stress, suggesting that the largest pores compacted the most.



**Figure 3:** The porosity of the pore networks of shale 4 as a function of pressure

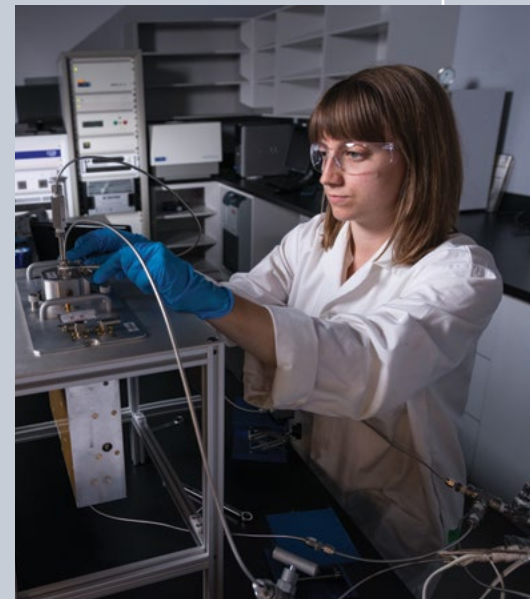

Finally, the right panel of Figure 2 shows the cumulative porosity of shale 2 as a function of confining pressure. Again it is obvious that the micro pores contribute to the majority of the total porosity of shale 2 while the macro pores contribute the least. Figure 3 shows the porosity data for shale 4 analyzed in the same fashion as the data for shale 2.




# Understanding changes in porosity as a function of overburden pressure in shales

Application Note 5



## Conclusion

A method has been presented for using NMR  $T_2$  relaxation times to identify pore networks in shales and their dependence on stress. Four samples were tested, all taken from wells in potential hydrocarbon-bearing formations. Two samples had a complex (trimodal) pore size distribution, and two contained only very small pores.

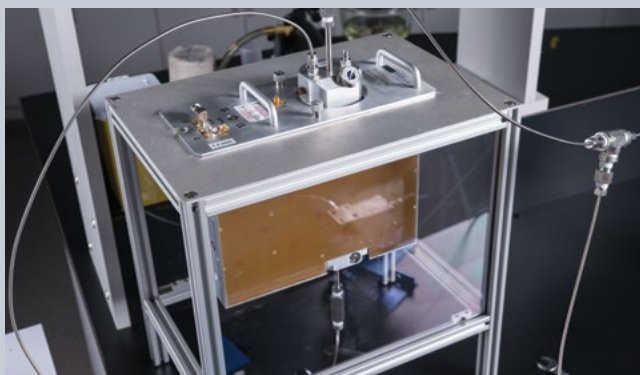

Where a complex distribution was present, the stress-dependence varied by pore size. The largest (macro) pores compressed significantly at relatively low stress (1,000 psi), and very little thereafter. This behavior is likely due to the closure of natural or induced fractures. The smallest (micro) pores show no significant change with stress. This is interpreted to indicate that these pores are not part of the mechanical fabric of the rock. The intermediate size (meso) pores show a continuing compression across the range of stresses.

## References

[1] Chhatre, S.S., Braun, E.M., Sinha, S., Determan, M.D., Passey, M.D., Zirkle, T.E., Wood, A.C., Boros, J.A., Berry, D.W., Leonardi, S.A. and Kudva, R.A., *Steady State Stress Dependent Permeability Measurements Of Tight Oil Bearing Rocks, 2014 Annual Symposium Of The Society Of Core Analysts (SCA2014-012)*.



## Oxford Instruments Magnetic Resonance

For more information: [magres@oxinst.com](mailto:magres@oxinst.com)  
[www.nmr.oxinst.com](http://www.nmr.oxinst.com)

### UK

Tubney Woods, Abingdon, Oxfordshire, OX13 5QX, UK  
**Tel:** +44 (0) 1865 393 200 **Fax:** +44 (0) 1865 393 333

### USA

300 Baker Avenue, Suite 150, Concord, MA, 01742, USA  
**Tel:** +1 978 369 9933 **Fax:** +1 978 369 8287

### China

Floor 1, Building 60, No.461, Hongcao Road,  
Shanghai, 200233, China  
**Tel:** +86 21 6073 2925 **Fax:** +86 21 6360 8535

## Green Imaging Technologies

For more information: [info@greenimaging.com](mailto:info@greenimaging.com)  
[www.greenimaging.com](http://www.greenimaging.com)

### Canada

520 Brookside Drive, Suite B, Fredericton, NB, E3A 8V2,  
Canada  
**Toll Free:** +1 888 944 8462  
**Tel:** +1 506 458 9992 **Fax:** +1 506 458 9615

This publication is the copyright of Oxford Instruments plc and provides outline information only, which (unless agreed by the company in writing) may not be used, applied or reproduced for any purpose or form part of any order or contract or regarded as the representation relating to the products or services concerned. Oxford Instruments' policy is one of continued improvement. The company reserves the right to alter, without notice the specification, design or conditions of supply of any product or service. Oxford Instruments acknowledges all trademarks and registrations.  
© Oxford Instruments plc, 2017. All rights reserved. Ref: MR/204/0717

