

# Benchtop NMR for Educational Experiments


## Nucleophilic Aromatic Substitution Reaction

This experiment is designed to teach the practical aspects and principles of the Nucleophilic Aromatic Substitution reaction:


### Refluxing, recrystallisation.

**Total experiment time: 2 hours**

A nucleophilic aromatic substitution is a reaction in which a nucleophile displaces a good leaving group, such as a halide, on an aromatic ring, and is of the general equation:

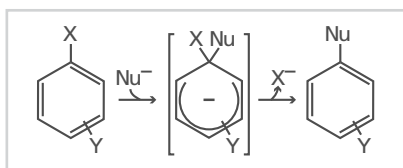

An example of the reaction is between 1-bromo-2,4-dinitrobenzene and aniline to form 2,4-dinitrodiphenylamine:

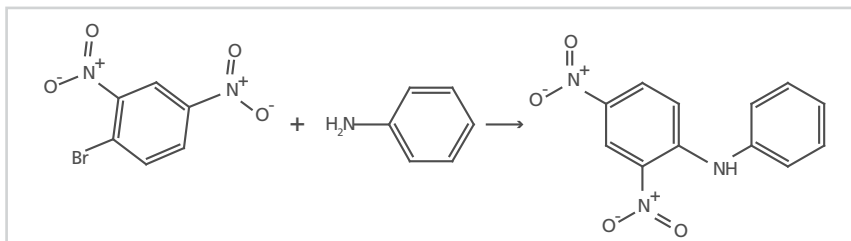

**Pulsar™**, a 60MHz benchtop NMR spectrometer can be used to measure the spectra of the starting materials and reaction product confirming that the reaction has completed and no starting materials remain. High quality spectra can be collected within 2 minutes. The spectra of the starting materials 1-bromo-2,4-dinitrobenzene and aniline, along with the final product 2,4-dinitrodiphenylamine are shown in Figure 1 (see over page).

It is clear from these spectra that 1-bromo-2,4-dinitrobenzene and aniline are not present in the spectrum of the product. The 2,4-dinitrodiphenylamine sample has tetramethylsilane (TMS) added as a chemical shift reference material showing a peak at 0ppm chemical shift.

The student can assign the peaks in the spectrum and generate peak integrals in order to verify the identity of the reaction product.

**OXFORD**  
INSTRUMENTS

*The Business of Science®*

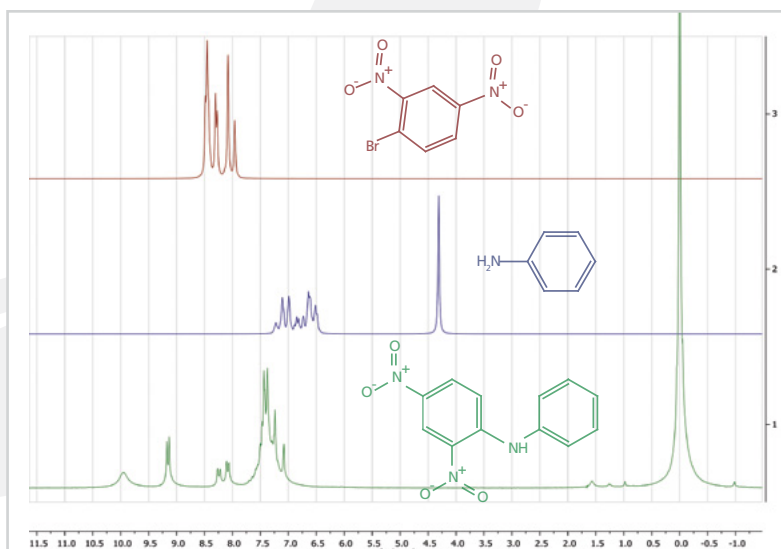


Figure 1. From top to bottom: 1-bromo-2,4-dinitrobenzene, aniline and 2,4-dinitrodiphenylamine spectra

**Pulsar** is a cryogen-free benchtop NMR spectrometer that can easily be sited in the undergraduate chemistry laboratory allowing hands-on NMR for students.

A simple to use software interface and standard sampling using 5mm NMR tubes enables a high throughput of samples in a busy laboratory.

Contact us now for more information

magres@oxinst.com



for more information visit [www.oxford-instruments.com/pulsar](http://www.oxford-instruments.com/pulsar)

This publication is the copyright of Oxford Instruments plc and provides outline information only, which (unless agreed by the company in writing) may not be used, applied or reproduced for any purpose or form part of any order or contract or regarded as the representation relating to the products or services concerned. Oxford Instruments' policy is one of continued improvement. The company reserves the right to alter, without notice the specification, design or conditions of supply of any product or service. Oxford Instruments acknowledges all trademarks and registrations. © Oxford Instruments plc, 2018. All rights reserved. Part no: Pulsar/Expt1/03-18



*The Business of Science®*